


PRODUCT INFORMATION


- 1** Flat lateral design
- 2** Smooth radius transition from proximal to distal portion of the nail
- 3** Flutes for stress reduced insertion
- 4** Distal diameter: 9, 10, 11 and 12 mm
- 5** CCD angle: 125° and 130°
- 6** Nail length: 170 mm, 200 mm, 240 mm left & right, 300 – 480 mm left & right



HERACLES PROXIMAL FEMUR NAIL

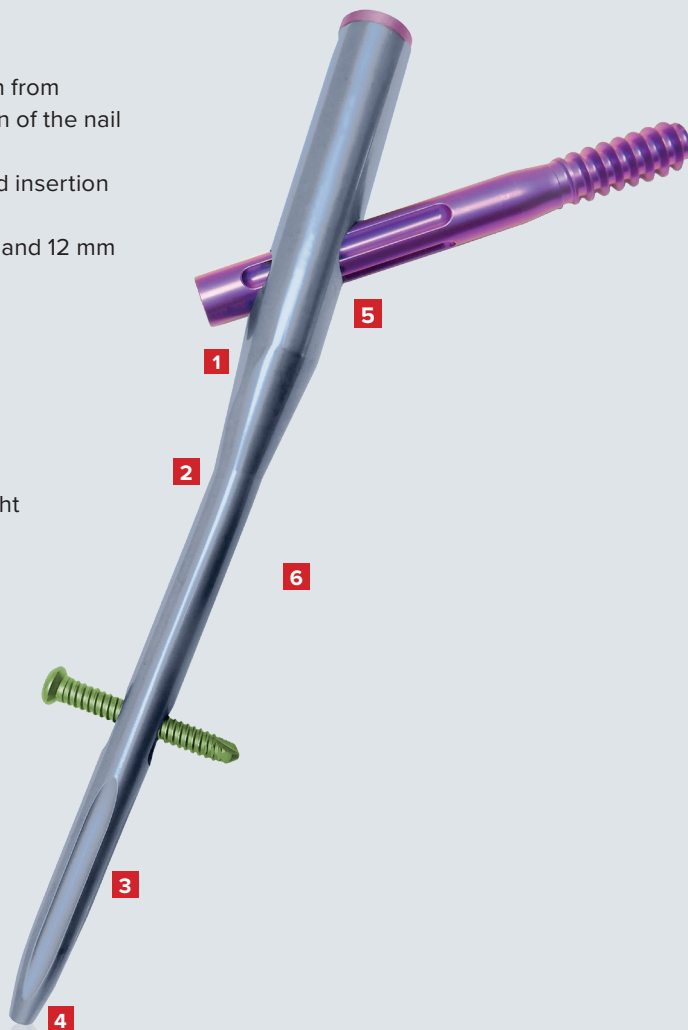

HERACLES PROXIMAL FEMUR NAIL

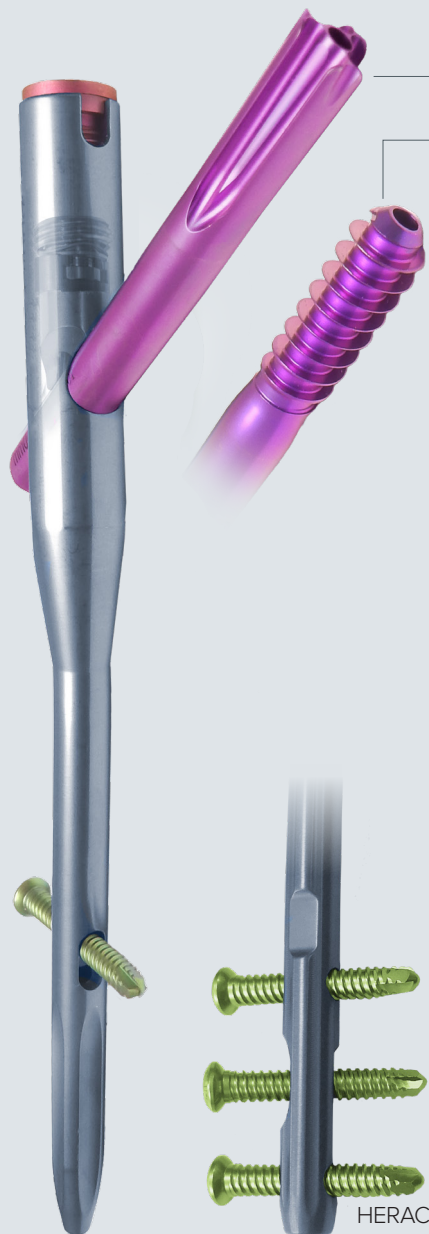

One System Providing Patient Specific Treatment Options

- Blade with three anchors for stable fixation in osteoporotic bone.
- Lag Screw with proven thread design.

Straightforward and Accurate Insertion


One piece radiolucent aiming arm:

- Less weight.
- Increased width to compensate excessive soft tissue.
- Aiming support possibilities.




HERACLES Long Nail

Smooth Nail Insertion to Preserve Reduction

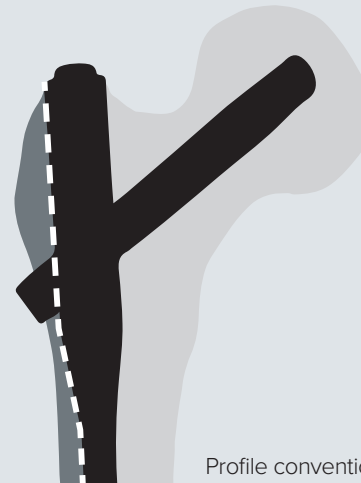


Profile conventional nail

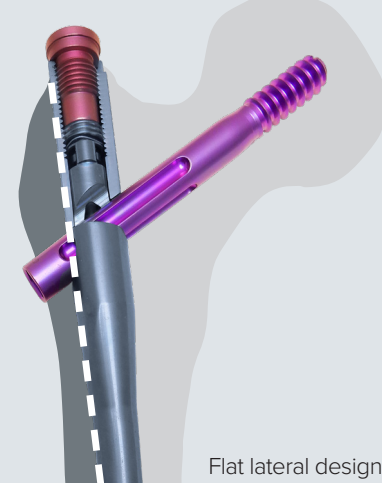


Flat lateral design

Just Lock with TwinLock

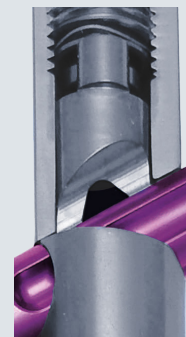

- Preassembled locking mechanism.
- Controlled and limited gliding.
- Locking of rotation.

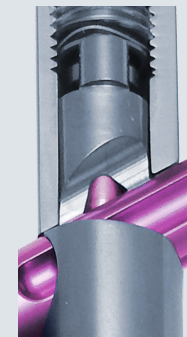

Unique End Cap for the option to stop gliding for certain indications.